



- Do **not** use tissue or toilet paper to clean the sensor, as it may scratch it.

- Not for use during or after swimming, bathing, or intense physical activity.

- Let the thermometer and person stay in the same environment (15°C–35°C / 59°F–95°F) for at least 30 minutes before measuring.

⚠ If problems persist, contact customer

- Expected use. The body temperature measured is displayed through the main display.




1. IR Sensor
2. Scan Button
3. Power-on Button
4. LCD Display
5. Battery Cover

- Automatically turn-off function (No operation 30 seconds)

1. Click “ ” to turn it on, the backlight will flash again.
2. 2 seconds later, it enter into the test state, meanwhile the blue backlight flashes.
3. Approach the thermometer to forehead at the distance of 0-3cm, then press “ ” to measure.



1. When the battery is low, the LCD will display the symbol, please replace the 3V CR2032 battery.
2. Remove the battery compartment cover by the screwdriver.

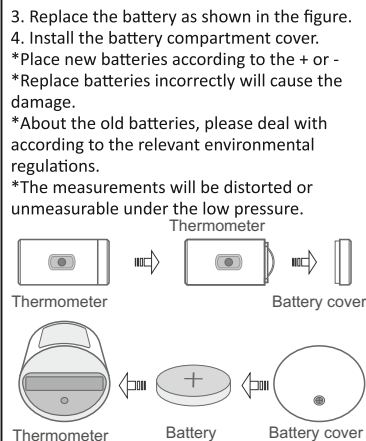

Register your product
to claim warranty



Register

* Terms & Conditions apply

 Manufactured by :
Passim Lifesciences Ltd.
Plot No. 45, Ind. Area, Phase-II, Panchkula - 134113, Haryana-INDIA
Mfg. Lic. No. : MFG/MD/2022/000615
For any complaint/Suggestion please contact :
Customercare Number : 1800 309, Timing : 9am - 7pm, Mon. - Sat.
Email ID : customercare@drodin.in, Website : www.drodin.in